

EXHIBIT A

[Login](#) [Register](#) [My Profile](#)
SIGMA-ALDRICH®
[Home](#) [Products](#) [Order Center](#) [Custom Products](#) [Support](#)

Product Name or No.

[Advanced Search](#)

[Product Catalog](#) > [Analytical/Chromatography](#) > [Resins and Media](#) > [Resins & LC Separation Media](#) > [Gel Filtration Media](#) > [Sephacore®](#)

Analytical / Chromatography

Products

> Air Monitoring

> Analytical Reagents

> Analytical Standards

> Gas Chromatography

> HPLC & LC/MS

> Microbiology

> Microscopy

> Sample Preparation & Purification

> Spectroscopy

> Syringes

> Titration

> Thin Layer Chromatography

> Vials & Accessories

> Videos

> Events - Trade Shows & Seminars

> Analytical Newsletters

> Reporter Europe 2008

Newsletter & e-News Registration

2008-09 Literature Request

> The Reporter

> Fluka - Analytical

> Supelco - Analytical

> Literature


 Sign-up for newsletters
& e-newsletters!

Analytix Reporter eAnalytix eReporter

Sign up today for the latest news, product information, and more!

Sephacore®

Product #	Description
25300	Sephacore® 2B 60-200 µm bead diameter
48200	Sephacore® 4B 40-165 µm bead diameter
68100	Sephacore® 6B 6% Beaded Agarose, 40-165 µm (wet), fractionation range 10,000-1,000,000 Da (dextrans)
CL2B300	Sephacore® CL 2B Cross-linked
CL4B200	Sephacore® CL 4B Cross-linked
CL6B200	Sephacore® CL 6B Cross-linked

[Reset Table](#)

[Site Use Terms](#) | [Terms and Conditions of Sale](#) | [Privacy](#) | [Contact Us](#) | [Site Map](#)

Copyright © 2008 Sigma-Aldrich Co. All Rights Reserved. Reproduction of any material from this site is strictly prohibited without permission. Sigma-Aldrich brand products are sold exclusively through Sigma-Aldrich, Inc.



SIGMA-ALDRICH


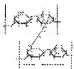
3050 Spruce Street
Saint Louis, Missouri 63103 USA
Telephone 800-325-5832 • (314) 771-5765
Fax (314) 286-7828
email: techserv@sial.com
sigma-aldrich.com

Product Information

SEPHAROSE AND SEPHAROSE CL GEL FILTRATION MEDIA

Exact replacement for Product Code 84963

The information below is abstracted, for customer convenience, from a supplier technical manual and a supplier information sheet. Complete copies of supplier data are available upon request from Sigma Research Technical Service.

SEPHAROSE	SEPHAROSE CL
 <p>Partial structure of Sepharose CL</p>	 <p>Partial structure of Sepharose CL</p>
CHARACTERISTICS	
<ul style="list-style-type: none">- Broad fractionation range- High exclusion limits- Negligible non-specific adsorption- Appearance: white suspension	<ul style="list-style-type: none">- Broad fractionation range- High exclusion limits- Negligible non-specific adsorption- Excellent chemical and physical stability- Appearance: white suspension
PRODUCT DESCRIPTION	
<p>Sephrose is a beaded agarose gel filtration medium with a broad fractionation range. Three different agarose contents are available: 2%, 4% and 6%, designated 2B, 4B and 6B, respectively. As agarose concentration increases porosity decreases, thus increasing rigidity and altering the fractionation range; nucleic acids and polysaccharides with molecular weights up to $\sim 4 \times 10^7$ can be separated on Sephrose 2B.</p>	<p>Sephrose CL is a cross-linked derivative of Sephrose, prepared by reacting Sephrose with 2,3-dibromopropanol under strongly alkaline conditions. After cross-linking, the gel is desulfated by alkaline hydrolysis under reducing conditions. The resulting crosslinked polysaccharide chains have a very low content of ionizable groups and better chemical and physical resistance than Sephrose, as well as improved flow properties. Porosity, however, is comparable to that of Sephrose. Three different agarose contents are available: 2%, 4% and 6%, designated CL-2B, CL-4B and CL-6B.</p>

				respectively.			
Property	2B	4B	6B	Property	CL-2B	CL-4B	CL-6B
MW Range:				MW Range:			
Globular Proteins	7×10^4 - 4×10^7	6×10^4 - 2×10^7	1×10^4 - 4×10^6	Globular Proteins	7×10^4 - 4×10^7	6×10^4 - 2×10^7	1×10^4 - 4×10^6
Dextrans	1×10^5 - 2×10^7	3×10^4 - 5×10^6	1×10^4 - 1×10^6	Dextrans	1×10^5 - 2×10^7	3×10^4 - 5×10^6	1×10^4 - 1×10^6
DNA exclusion limit	-1352 base pairs	-872 base pairs	-45-165 base pairs	DNA exclusion limit	-1353 base pairs	-872 base pairs	-45-165 base pairs
Diameter (wet bead)	-60-200 μ m	-45-165 μ m	-45-154 μ m	Diameter (wet bead)	-60-200 μ m	-45-165 μ m	-45-165 μ m
pH Range	4-9	4-9	4-9	pH Range	3-14	3-14	3-14
Max. Pressure	40 cm H ₂ O	80 cm H ₂ O	200 cm H ₂ O	Max. Pressure	50 cm H ₂ O	120 cm H ₂ O	>200 cm H ₂ O
Max. Volumetric Flow Rate	0.83 mL/min	0.96 mL/min	1.16 mL/min	Max. Volumetric Flow Rate	1.25 mL/min	2.17 mL/min	2.5 mL/min
Max. Linear Flow Rate	10 mL/cm ² h	11.5 mL/cm ² h	14 mL/cm ² h	Max. Linear Flow Rate	15 mL/cm ² h	26 mL/cm ² h	30 mL/cm ² h

STABILITY

Sephacrose melts upon heating to 40°C, cannot be autoclaved, and the bead structure may be damaged upon freezing. Due to the presence of 3,6-anhydro-L-galactose, the matrix is resistant to biological degradation. Sephacrose is stable in aqueous (including saline) solutions at pH 4-9. Use of dissociation media such as guanidine hydrochloride and urea, chaotropic salts such as KSCN, and oxidizing agents is not advisable because these reagents may disrupt the hydrogen bonds which stabilize the matrix. Sterilization can be done chemically (e.g., by treatment with DEPC).

Sephacrose CL is stable in aqueous media at pH 3-14, and in organic solvents. Dissociating media and chaotropic salts may be used, but oxidizing conditions should be avoided in order to preserve the identity of the constituent sugar residues. Sephacrose CL may be autoclaved repeatedly at 120°C and pH 7. The matrix is resistant to biological degradation. Flow rates may be up to 50% higher than on Sephacrose with an equivalent fractionation range.

USAGE AND REGENERATION

Sephacrose and Sephacrose CL are supplied pre-swollen as suspensions in distilled water. Before packing a column, dilute the required amount of gel with starting buffer to form a thick slurry, about 75% of which is settled resin, then degas the slurry. The maximum operating pressures (shown above) should never be exceeded. Pass 2-3 column volumes (CV) of eluent through the gel to equilibrate the bed. Elution may be achieved by gravity feed or through the use of a peristaltic pump; more consistent flow rates and more reproducible separations are obtained with a pump. Sephacrose contains a small number of ionic sulfate and carboxyl groups which may cause adsorption of basic proteins at low ionic strengths. Therefore, eluents with

ionic strengths exceeding 0.02 M are sometimes necessary (for example, tRNA species have been resolved in high concentrations of ammonium sulfate, and DNA has been separated from RNA in 1.5 N NaCl). The gels can be cleaned as indicated below and stored at 4-8°C in a suitable antimicrobial agent (e.g., 20% ethanol) for indefinite time periods.

Sepharose should be cleaned in the column or batchwise with a non-ionic detergent solution.

Sepharose CL can be cleaned in the column or batchwise. Wash with at least 1 CV of 0.5 N NaCl in 0.1 N NaOH (a step which also sterilizes the resin), then with 10 CV of water (or until the eluent is at neutral pH).